



地址：上海市张江高科技园区爱迪生路326号 邮编：201203

销售热线：400-820-3059

86-21-51320588-8008/8035/8026/8025


传真：86-21-51320502

QQ:800015916

MSN:Tauto_shanghai@hotmail.com

URL: www.tautobiotech.com

Mail: sales2@tautobiotech.com export2@tautobiotech.com

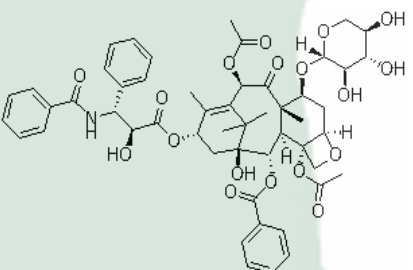PRODUCT ID	E-0520	Structure 
PRODUCT	7-xylosyltaxol/7-木糖基紫杉醇	
STANDARD TYPE	Quantitative Reference Material (含量测定)	
LOT NUMBER	18030124	
ASSAY METHOD	HPLC,MS,NMR	
DATE OF SAMPLE	2018-03-01	
DATE OF REPORT	2018-03-01	

OTHER NAME Benzenepropanoic acid, 7-(β -Xylosyl)taxo

CHEMICAL FORMULA C₅₂H₅₉NO₁₈

MOLECULAR WEIGHT(MW) 986.03

CAS NUMBER 90332-66-4

CHEMICAL FAMILY NA

ANALYTICAL RESULT

TEST	METHOD	SPECIFICATION	RESULT
HPLC	Q/TFFK319—2007	NLT $\geq$ 97.0%	97.3%
NMR		Conforms	Conforms
MS		Conforms	Conforms
Residual Solvent	Q/TFFK319—2007	NA	ND
Water	Q/TFFK319—2007	$\leq$ 2%	Conforms
Appearance		White Powder	Conforms

STORAGE CONDITIONS

STORAGE Refer to the label

SHELF LIFE 2 Years


Analyst

Checker

Manager

Product Analysis Report

Product Name: 7-xylosyltaxol **Chrom-Column:** C18
Lot Number: 18030124 **Mobile Phase:** ACN:H₂O=45:55
Temperature: 25°C **Wavelength:** 230nm
Detect Time: 2018-03-01 **Inject Vol:** 20µL
Sample: 70%MeOH **Flow:** 1ml/min

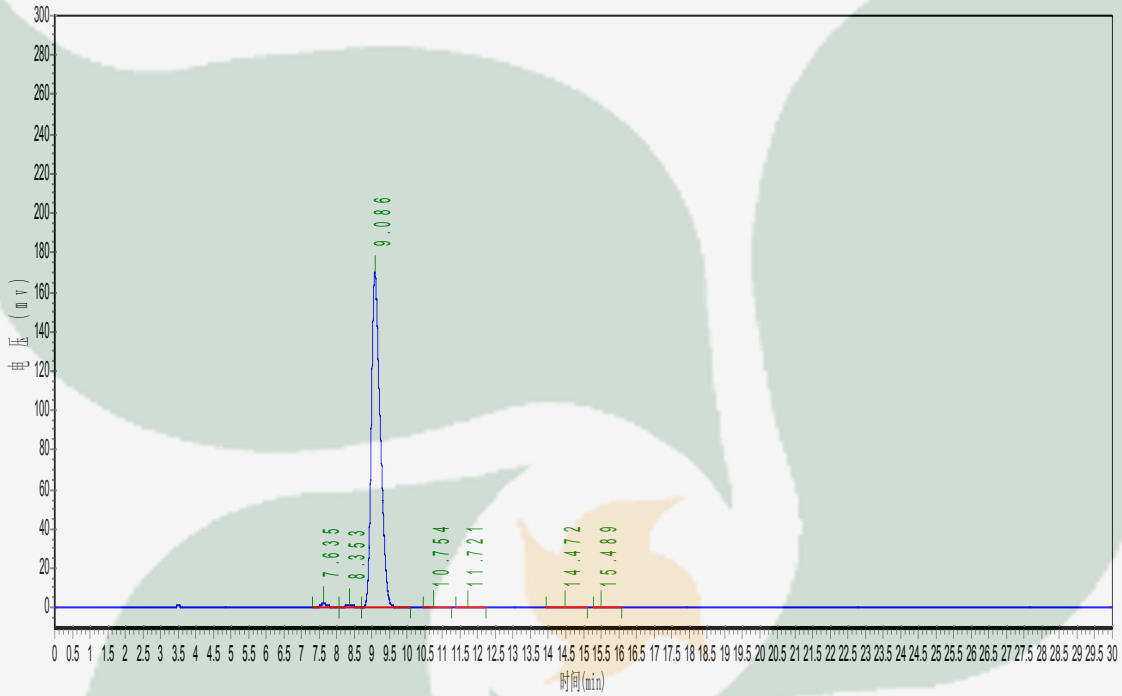

No.	Ret.Time min	Height mAU	Area mAU*min	Rel.Area %
1	7.635	2105.886	33100.520	1.0795
2	8.353	1599.545	25091.137	0.8183
3	9.086	170008.547	2983766.750	97.3059
4	10.754	342.102	6289.900	0.2051
5	11.721	375.944	8113.264	0.2646
6	14.472	149.514	4465.700	0.1456
7	15.489	216.102	5550.975	0.1810