



地址：上海市张江高科技园区爱迪生路326号 邮编：201203

销售热线：400-820-3059

86-21-51320588-8008/8035/8026/8025


传真：86-21-51320502

QQ:800015916

MSN:Tauto_shanghai@hotmail.com

URL: www.tautobiotech.com

Mail: sales2@tautobiotech.com export2@tautobiotech.com

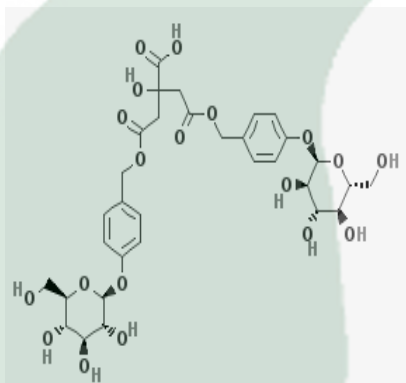PRODUCT ID	E-0817	Structure 
PRODUCT	Parishin C/巴利森苷C	
STANDARD TYPE	Quantitative Reference Material (含量测定)	
LOT NUMBER	18031421	
ASSAY METHOD	HPLC,MS,NMR	
DATE OF SAMPLE	2018-03-14	
DATE OF REPORT	2018-03-14	

OTHER NAME	柠檬酸酯 C
CHEMICAL FORMULA	C ₃₂ H ₄₀ O ₁₉
MOLECULAR WEIGHT(MW)	728.65
CAS NUMBER	174972-80-6
CHEMICAL FAMILY	Miscellaneous

ANALYTICAL RESULT

TEST	METHOD	SPECIFICATION	RESULT
HPLC	Q/TFFK542—2007	NLT ≥98.0%	98.1%
NMR		Conforms	Conforms
MS		Conforms	Conforms
Residual Solvent	Q/TFFK542—2007	NA	ND
Water	Q/TFFK542—2007	≤2%	Conforms
Appearance		white powder	Conforms

STORAGE CONDITIONS

STORAGE

Refer to the label

SHELF LIFE

2 Years

Analyst

Checker

Manager

Product Analysis Report

Product Name: Parishin C

Chrom-Column: C18

Lot Number: 18031421

Mobile Phase: ACN:0.1%H₃PO₄=15:85

Detector: UV


Wavelength: 220nm

Detect Time: 2018-03-14

Inject Vol: 20ul

Sample dissolution: H₂O

Flow: 1.0ml/min

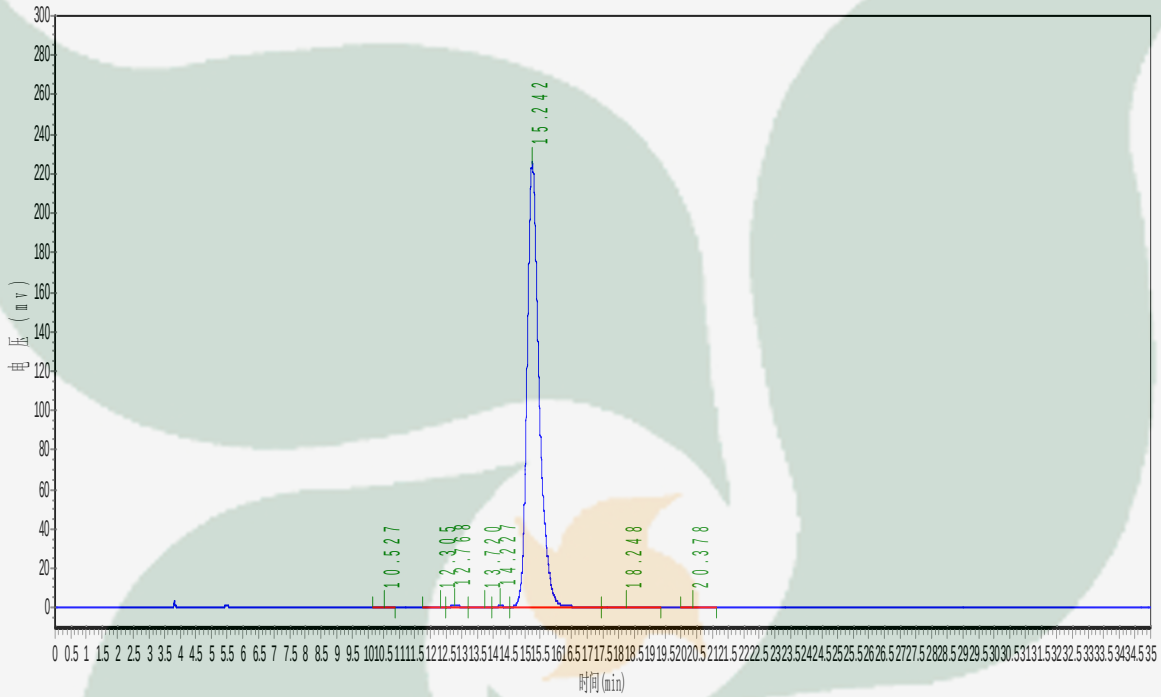

No.	Ret.Time min	Height mAU	Area mAU*min	Rel.Area %
1	10.527	340.239	4733.500	0.0762
2	12.305	394.878	7733.275	0.1246
3	12.768	1339.400	32478.025	0.5231
4	13.720	613.522	20005.523	0.3222
5	14.227	952.187	24402.514	0.3930
6	15.242	225620.516	6091651.500	98.1117
7	18.248	588.387	23646.418	0.3808
8	20.378	168.339	4240.150	0.0683