



地址：上海市张江高科技园区爱迪生路326号 邮编：201203

销售热线：400-820-3059


86-21-51320588-8008/8035/8026/8025

传真：86-21-51320502

QQ:800015916

URL: www.tautobiotech.com

Mail: sales2@tautobiotech.com export2@tautobiotech.com

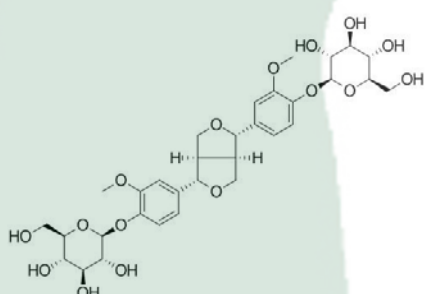PRODUCT ID	E-0227	Structure 
PRODUCT	Pinoresinol Diglucoside /松脂醇二葡萄糖苷	
STANDARD TYPE	Quantitative Reference Material (含量测定)	
LOT NUMBER	18021331	
ASSAY METHOD	HPLC,MS,NMR	
DATE OF SAMPLE	2018-02-13	
DATE OF REPORT	2018-02-13	

OTHER NAME	NA
CHEMICAL FORMULA	C32H42O16
MOLECULAR WEIGHT(MW)	682.66628
CAS NUMBER	63902-38-5
CHEMICAL FAMILY	Lignans

ANALYTICAL RESULT

TEST	METHOD	SPECIFICATION	RESULT
HPLC	Q/TFFK228—2007	NLT $\geq$ 97.0%	97.5%
NMR		Conforms	Conforms
MS		Conforms	Conforms
Residual Solvent	Q/TFFK228—2007	NA	ND
Water	Q/TFFK228—2007	$\leq$ 2%	Conforms
Appearance		White Powder	Conforms

STORAGE CONDITIONS

STORAGE	Refer to the label
SHELF LIFE	2 Years


Analyst

Checker

Manager

Product Analysis Report

Product Name: Pinoresinol Diglucoside Chrom-Column: C18
Lot Number: 18021331 Mobile Phase: MeOH:H₂O=30:70
Detector: UV Wavelength: 277nm
Detect Time: 2018-02-13 Inject Vol: 20ul
Analyst: shuzhilin Flow: 1.0ml/min

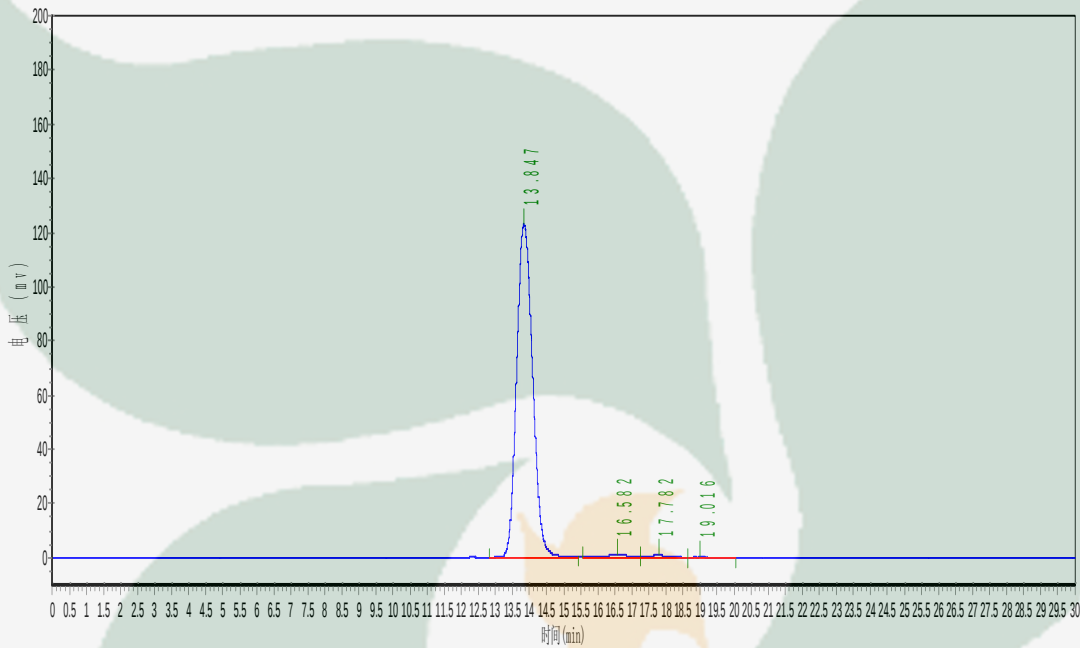

No.	Ret.Time min	Height mAU	Area mAU*min	Rel.Area %
1	13.847	122839.367	4018852.000	97.5082
2	16.582	1160.002	55785.363	1.3535
3	17.782	886.963	39864.570	0.9672
4	19.016	166.395	7050.224	0.1711